




**The professional concept for
home usage to lose weight**

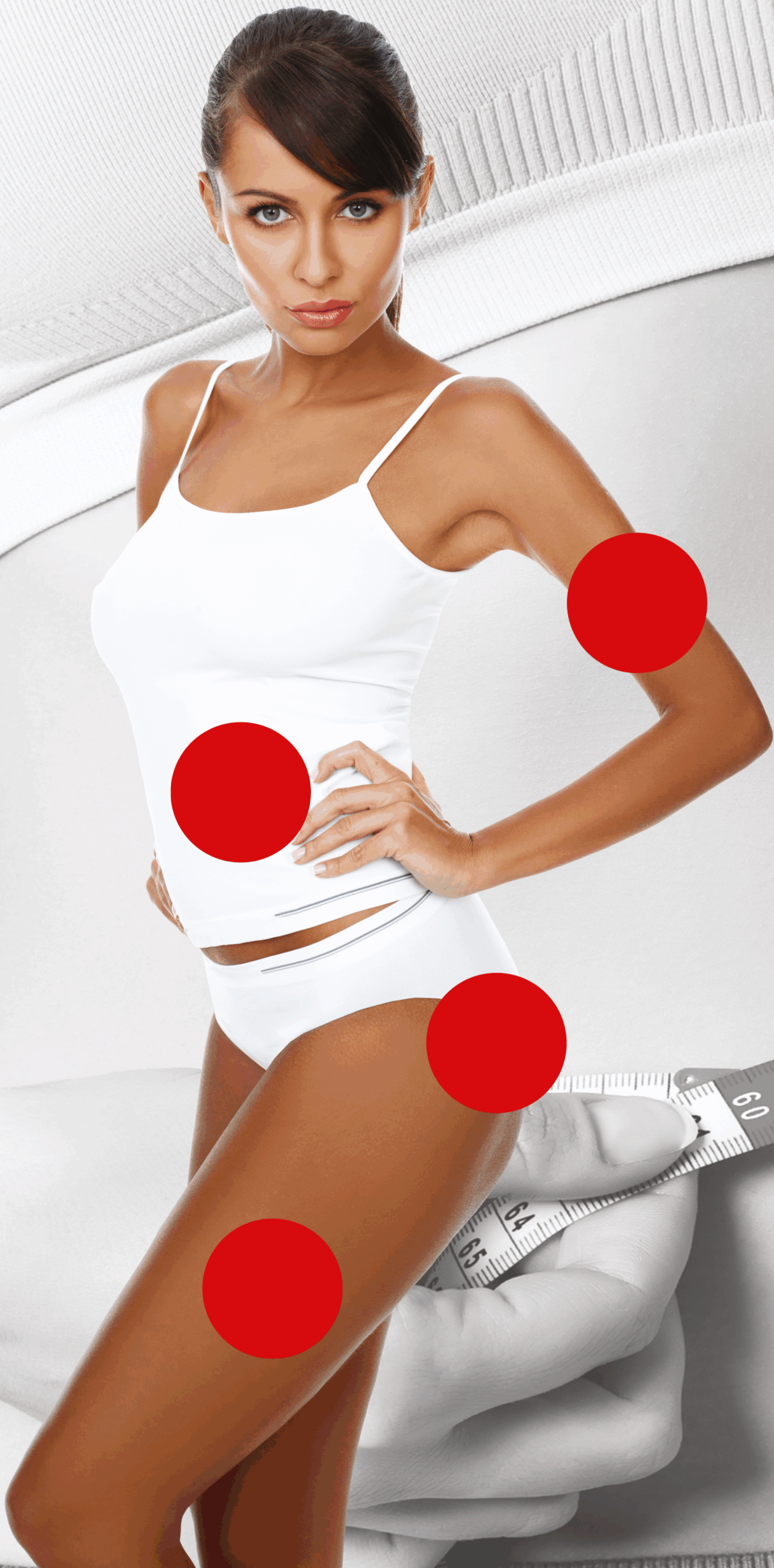

TRISLIM is a newly developed product with the goal to improve your outer appearance by pointing at the main problem areas:

belly - legs – arms – buttocks.

Our goal was to develop a system which:

- achieves maximum results
- causes few costs
- is easy to handle

The problem

Today it is assumed that so-called problem areas develop due to metabolic disorders caused by a change of life style, illnesses or lack of exercise.

All problem areas have three main factors in common:

1. a bad blood flow connected to a metabolic disorder
2. an increased accumulation of liquids in the connective tissue
3. fat cells with increased fat quantity

Therefore, an effective therapy must:

- improve the blood flow
- decrease the amount of liquids in the connective and fat tissue
- encourage the metabolism
- accelerate fat reduction



The solution TRISLIM

Working with air-controlled cuffs the system can massage a specific, large zone of problem areas.

The TRIMSLIM one-chamber system is directional, i.e. the compressed air and consequently also the massage are performed in defined directions which guarantee not only an optimization of the blood flow but also the best activation of the lymphatic system.

The working method is simple:



Taking the leg cuff as an example, venous blood and lymphatic fluid is transported to the body centre during the pressure phase. The reduction of fat is increased in this region due to an increased blood flow. During the normal pressure phase arterial blood is increasingly transported to the capillaries.

Effect

The effect of TRISLIM is similar to the one of manual lymph drainage only that this one is performed on a bigger area and for a longer time. TRISLIM works with fast changing overpressure. The cuffs are inflated in a controlled way thus applying pressure on the connective/ fat tissue. During the pressure phase venous blood and lymphatic liquid is transported away. During the next phase with normal pressure arterial blood can flow. Contrary to similar simple devices, using TRISLIM you can set the time and pressure intervals. Thus, you can individually adjust the system to the tissue structure of each customer.

Important for a good effect is:

1. the pressure set
2. the time interval set

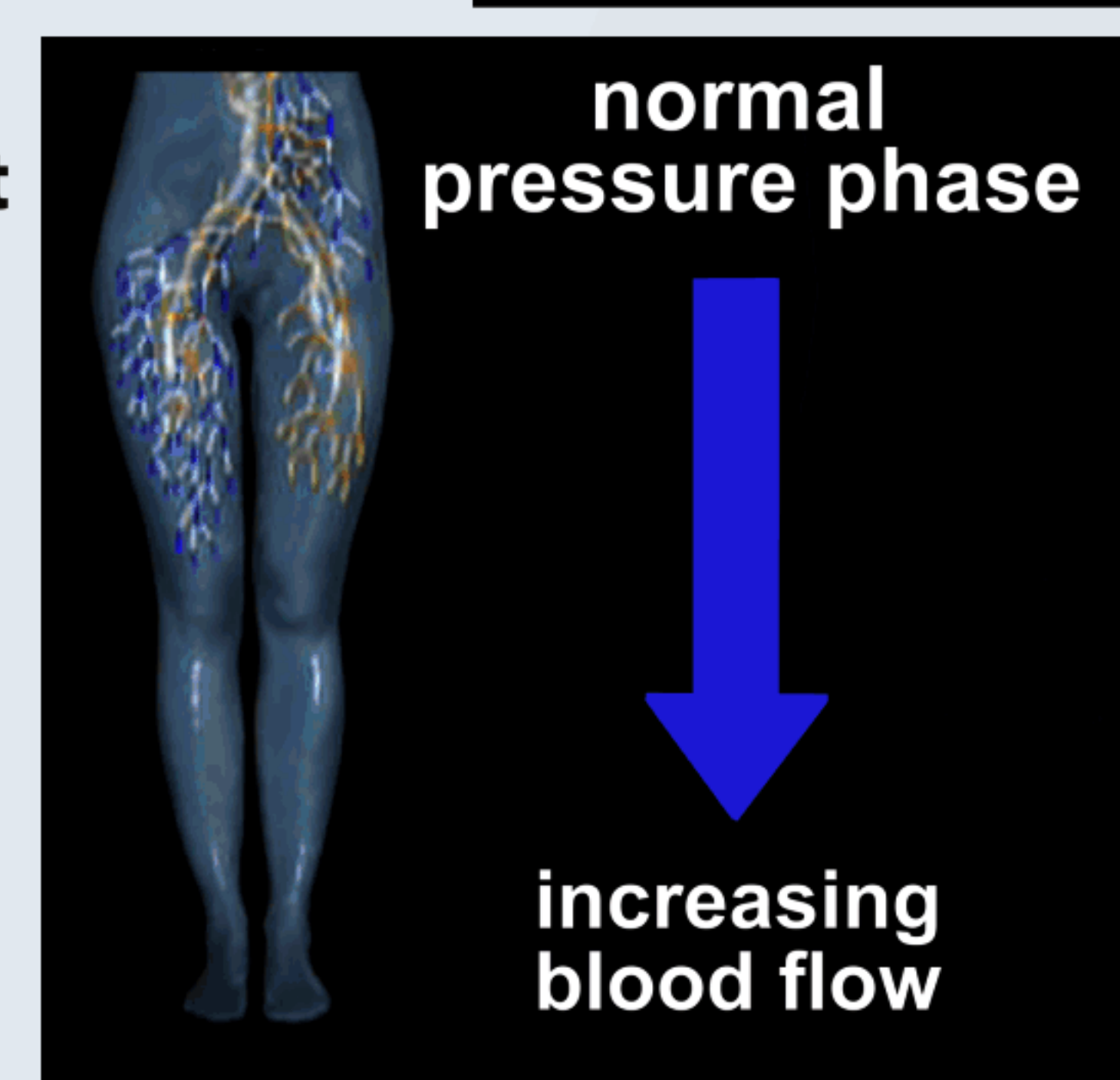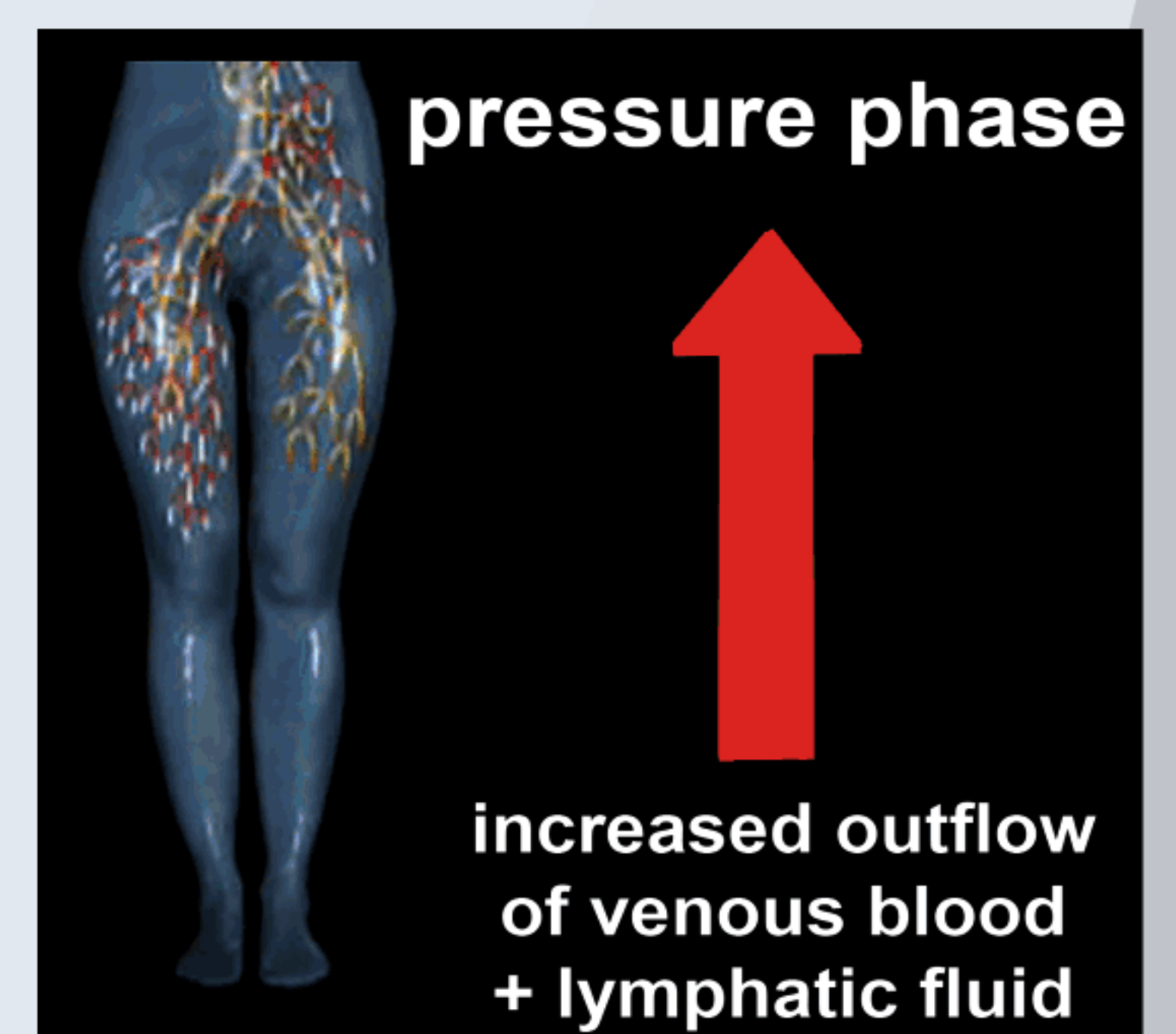


Practice

TRISLIM is perfect in connection with normal cardio training.
 The increased blood flow, the "muscle-pumping effect" and the massage of TRISLIM complete each other during local size reduction.
 One training unit with TRISLIM is 30 minutes.
 Basically every cardio device is suitable for a training with TRISLIM.
 The system can be handled easily and the cuffs can be put on by your customers.
 The chip card contains the credit of training units.

What are the results ?



Many physiological values must be considered in order to make exact prognoses about a customer. However, these values may be available in rare cases only.
 A size reduction of 3-7 cm in 4 weeks is realistic given a certain disposition and 3 training units a week. Simultaneous nutrition counselling together with light weight training increase the training effect.
 Before/ after 6 weeks = 18 units

Result - 6,0 / -5,7 cm | r / l

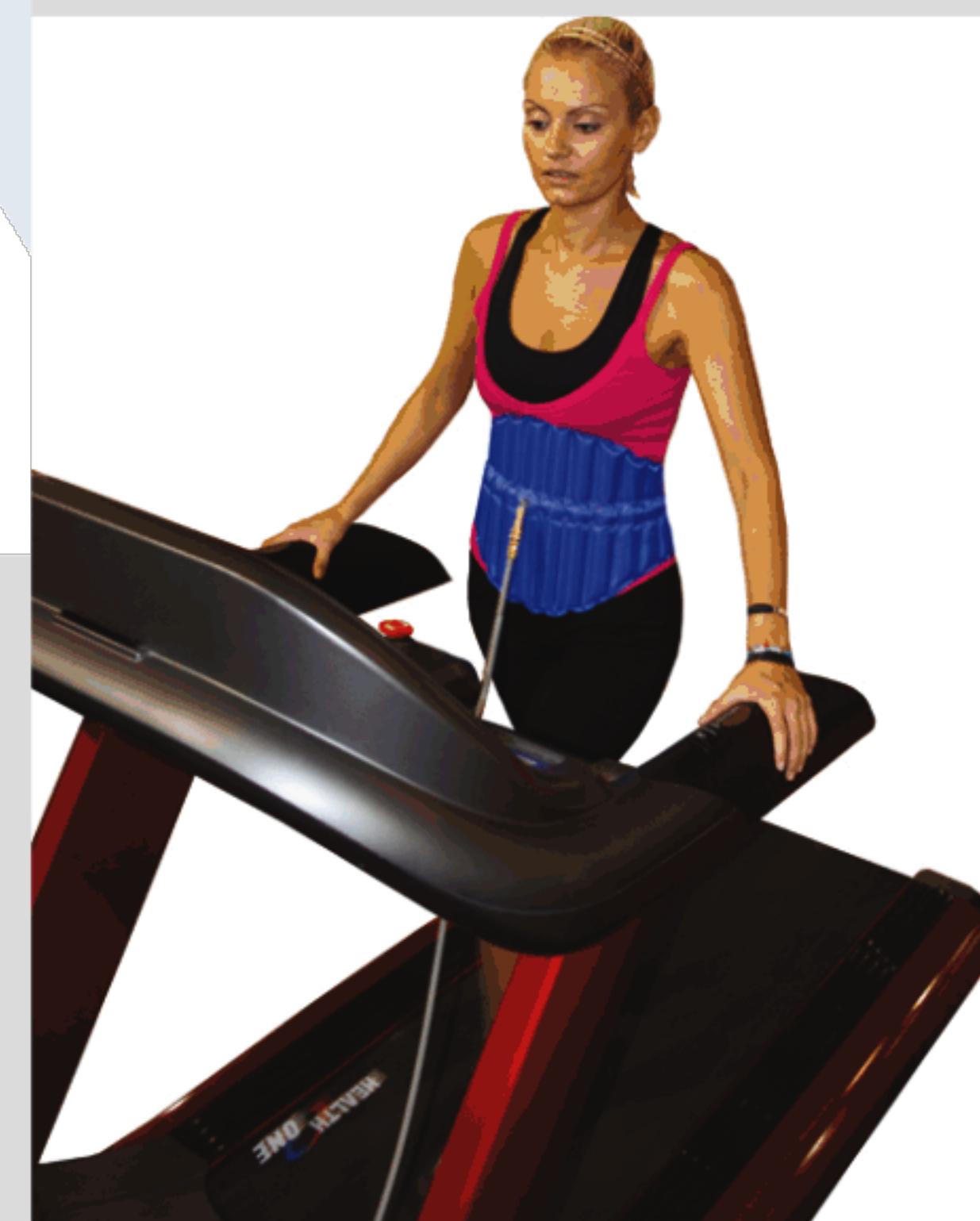

Cuffs

TRISLIM cuffs are made of two different materials. The outer blue material is made of PU plastic which is extremely resistant against breaking and rubbing.
 The inner part is made of pure PU: a very stretchy material that easily adapts to the body.
 Thanks to this material the different pressure levels can be perfectly transferred to the body. PU is skin-friendly, easily to disinfect and also often used in medical technology. We have chosen this material combination in order to guarantee not only a big lifetime but also a good adaption to the body. The cuffs are made in Germany in a factory certified according to DIN 9001 and ISO EN 13485.

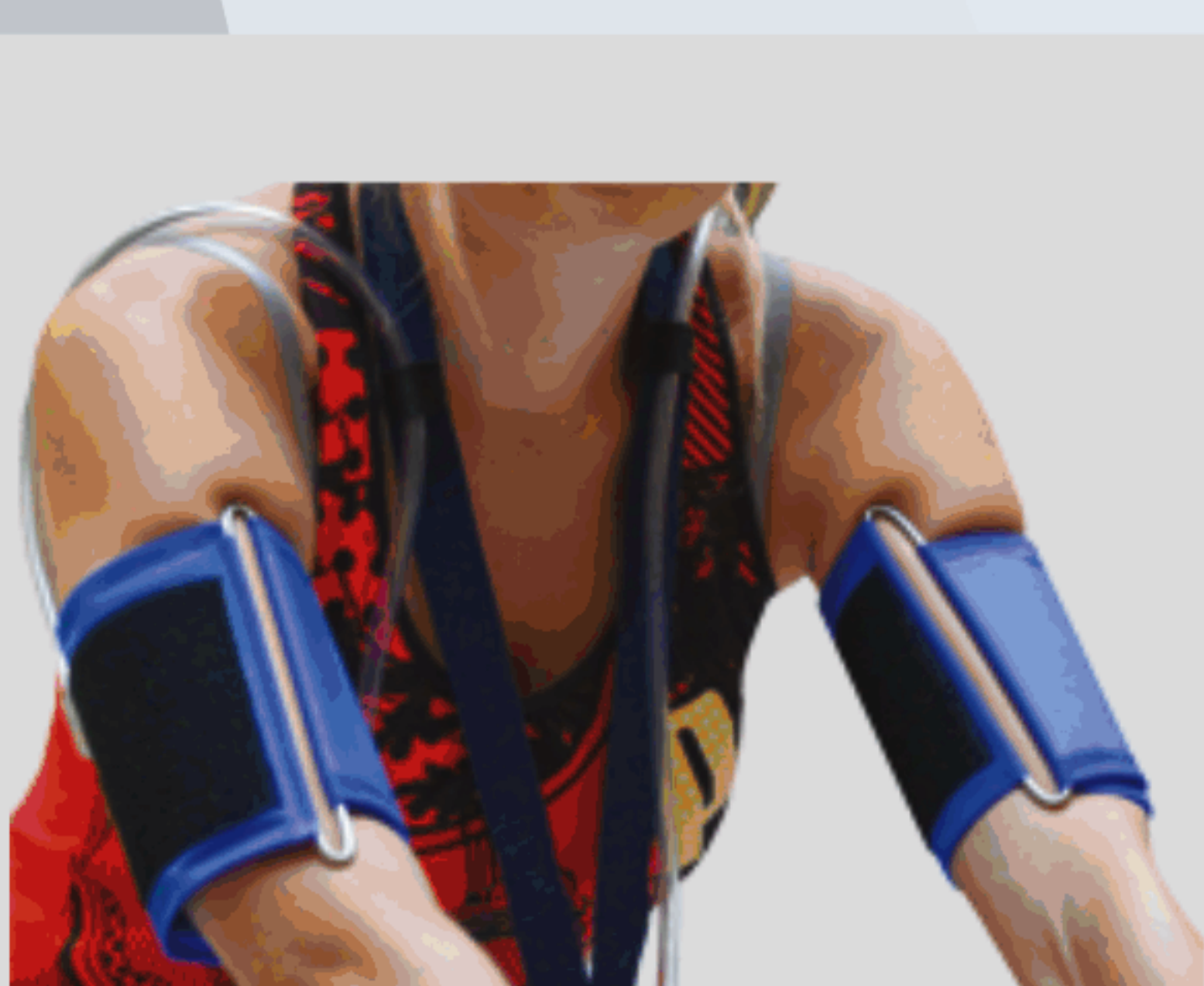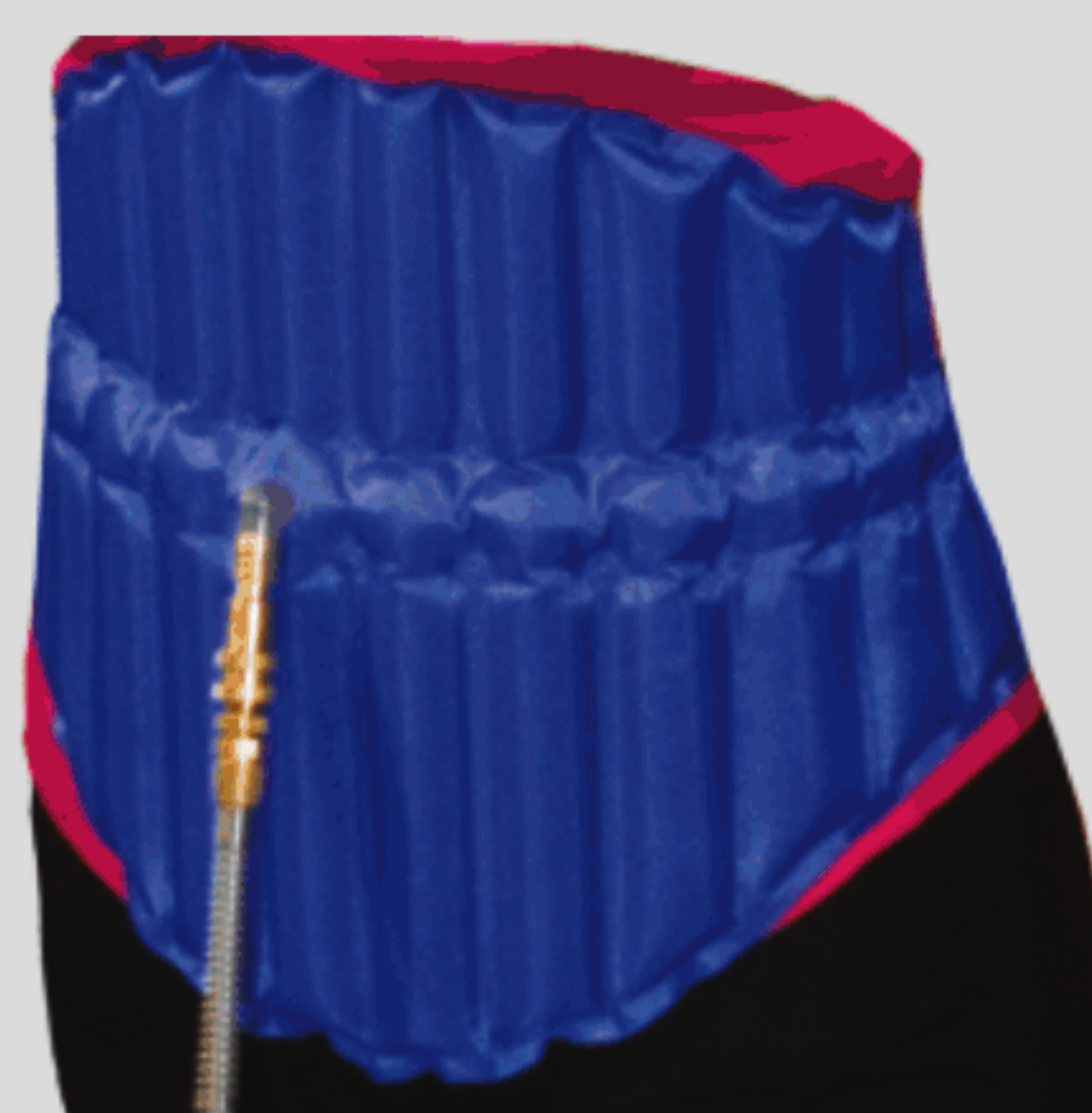



4 Problem areas 1 Solution



Technical specifications:

Control Unit :

pressure:
time intervall:
dimensions (w x h x d) :
weight :
power connection
with separate AC adapter:
connection hose :

LCD -Display 240 x 128 Px
2- key control
microprocessor controlled

adjustable to 5 levels, 150-200-300-400-500 mbar
20 sec.
220 x 950 x 210 mm
3,00 Kg
100 - 240 V, ~ 1,5 A, 50-60HZ
2,50 m, silicone, connection over fast coupling



Cuffs

material: outside: resistant nylon textile covered with PU
inside - PU

quick release fastener with plug fastener, belt PE 40 mm

Belly cuff : size 1 for 55 - 95 cm(girth)
size 2 for 85 - 145 cm(girth)

Leg cuff: for 55 - 85 cm (girth)
Delivery in pairs (measured at the biggest point)

Arm cuff : one size 20 - 40 cm
Delivery in pairs

Buttocks cuff: one size 90 - 140 cm

TRISLIM is made in Germany